



12

FASTRAK

ALL PRODUCTS SHIP IN
FIVE BUSINESS DAYS

SCIENCE | LAB | RESEARCH | TECH ED | ART | INDUSTRIAL ARTS



TABLES



PLAIN APRON TABLES—Legs and aprons are constructed of solid red oak, and are finished with a chemical resistant, earth-friendly UV finish. Aprons are made of solid red oak and are 13/16"D x 4-1/2"H with heavy-duty 13 gauge steel corner braces. 2-1/4" square solid oak legs, with 5/16" hanger bolts installed, are attached to apron with nuts and washers supplied. All tables are equipped with non-skid adjustable glides for precise leveling and for prevention of rocking on uneven floors. Tables are supplied with black rubber leg boots which protect against scuffing and chemical spills. Standard table height is 30", with a weight capacity of 500 lbs. (Shipped KD with tops attached; FOB Wisconsin, 54174)

SIZE	P. LAM	CHEMGUARD	SOLID EPOXY
24" x 48"	P7101K30N*	P7102K30N*	P7106K30N*
24" x 54"	P7201K30N*	P7202K30N*	P7206K30N*
24" x 60"	P7601K30N*	P7602K30N*	P7606K30N*
24" x 72"	P7301K30N*	P7302K30N*	P7306K30N*
30" x 60"	P7141K30N*	P7142K30N*	P7146K30N*
30" x 72"	--	P7152K30N*	--

MAPLE SCIENCE TABLES—Constructed of solid maple and finished with an earth-friendly UV finish. Aprons are constructed of solid maple and are 13/16"D x 4-1/2"H with heavy-duty 13-gauge steel corner braces. It has solid 2-1/4" maple legs with 5/16" hanger bolt installed, and are attached to the table with nuts and washers (supplied). All tables include leg boots and levelers. Standard height is 30"H. Shipped KD with tops attached.



PLAIN APRON	1-1/4" P. LAM	1" POWDER COAT	1-1/4" CHEMGUARD	1" SOLID EPOXY
48"W x 24"D	P7101M30N*	P7109M30N*	P7102M30N*	P7106M30N*
54"W x 24"D	P7201M30N*	P7209M30N*	P7202M30N*	P7206M30N*
60"W x 24"D	P7601M30N*	P7609M30N*	P7602M30N*	P7606M30N*

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



TABLES

BOOK COMPARTMENT TABLES—

Are built with the same quality techniques as the plain apron table. They are finished with a chemical resistant, earth-friendly UV finish. These tables have two book compartments installed in the apron. The book compartments measure 14-1/4"W x 12-3/4"D x 3-7/8"H. The table has a weight capacity of 500 lbs. (Shipped KD with tops attached; tables ship FOB Wisconsin, 54174)



SIZE	COMPARTMENTS	P. LAM	CHEMGUARD	SOLID EPOXY
24" x 48"	2	C7101K30N*	C7102K30N*	C7106K30N*
24" x 54"	2	C7201K30N*	C7202K30N*	C7206K30N*
24" x 60"	2	C7601K30N*	C7602K30N*	C7606K30N*
24" x 72"	2	C7301K30N*	C7302K30N*	C7306K30N*
30" x 60"	2	C7141K30N*	C7142K30N*	C7146K30N*
30" x 72"	2	--	C7152K30N*	--



ADJUSTABLE HEIGHT TABLE—A perfect solution for ADA requirements with the ability to adjust in 2 cm increments from 27"H to 39"H. Can also be adjusted for use as an incline surface for physics. Solid 4-1/2" oak aprons are attached to your choice of two top options. The wood is finished with a chemical resistant, earth-friendly UV finish. Optional casters sold separately. (FOB Wisconsin, 54174)



SIZE	P. LAM	CHEMGUARD
24" x 48"	P8101K*	P8102K*
24" x 54"	P8201K*	P8202K*
24" x 60"	P8601K*	P8602K*
24" x 72"	P8301K*	P8302K*
30" x 60"	P8141K*	P8142K*
Casters Option	100709X4* (Set of Four)	

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

WORKSTATIONS



FORWARD VISION I—The unique design leaves enough room to store four stools (not included) under the workstation. Sink units comes with two combination cold water/gas fixtures. It comes with two GFI electrical outlets, and ample storage. The stations are finished with a chemical resistant, earth-friendly UV finish. All doors and drawers lock. Dimensions: 96"W x 50"D x 36"H. (FOB Wisconsin, 54174)

ITEM #	SINK DIMENSIONS	TOP
2946K*	12"W x 8"D x 6"H	Solid Epoxy

FORWARD VISION II—Teachers maintain control of the class with all students facing the instructor. The large work surface—46"D x 88"W—is made of chemical resistant phenolic resin and holds the solid epoxy sink, and two combination cold water & gas fixtures. The cabinets provide locked storage space and side mounted hooks hold book bags off the floor. It is finished with a chemical resistant, earth-friendly UV finish. Includes two GFI AC duplex units. This model is compact enough to fit multiple units in a small classroom and is easy to assemble. Dimensions 88"W x 46"D x 36"H. Wt. 503 lbs. (FOB Wisconsin, 54174)



ITEM #	SINK DIMENSIONS	TOP
2944K*	12"W x 8"D x 6"H	Solid Phenolic

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



WORKSTATIONS

CLEAN UP SINK—One piece top includes sink of molded polyolefin (easy to clean, extremely chemical resistant, and is highly impact resistant to reduce glassware breakage). Long grooved slopping counter leads to large sink (24-1/2"W x 18-1/2"D x 14"H) in center. Unit includes a hot and cold water mixing faucet. Constructed of hardwood and oak veneers, with a chemical resistant, earth-friendly UV finish. Two adjustable shelves can be found inside the outer cabinets. It includes pulls, locks and rubber base molding. (FOB Wisconsin, 54174)



ITEM #	SIZE	TOP
3303K*	55-1/2"W x 28"D x 36-1/2"H	Polyolefin



MIRROR—The shatter resistant acrylic mirror is enclosed in a solid oak frame with a natural transparent finish. Easily mounts on 1/2" diameter crossbar of any of our mobile labs which include uprights and crossbar. It can be tilted to any angle—held with nylon screws in the mounting brackets. Weight 10 lbs. (FOB Wisconsin, 54174)



MIRROR/MARKER BOARD—This solid oak frame unit allows for the unobstructed use of the shatter resistant acrylic mirror and marker board. Unit can be tilted to any angle or completely reversed while the demonstration is in progress. The unit is held secure to the 1/2" diameter crossbar by means of screws in the frame. Weight 12 lbs. (FOB Wisconsin, 54174)

ITEM #	SIZE	PRODUCT
4000K*	22-1/2" x 34-1/2"	Mirror
4001K*	20-3/4" x 27-3/4"	Mirror/Marker Board

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

INSTRUCTOR'S DESKS



5' INSTRUCTOR'S DESK—Constructed of solid oak, oak hardwood and oak veneers, this unit serves as a complete demonstration center. It is finished with a chemical resistant, earth-friendly UV finish. Included are four drawers (three small, one file), and a large cupboard. The drawers are expertly crafted from the finest materials with dovetail joinery throughout. Drawers run on epoxy-coated steel glides with nylon bearings for smooth operation. Includes a combination cold water/gas fixture and a epoxy sink that measures 14"W x 10"D x 6"H. All units come with a GFI protected AC duplex receptacle. Dimensions; Knee space: 22"W x 32-1/2"D; Pull-out: 3/4"H x 16-1/4"W x 25-1/4"D; Overall: 60"L x 30"D x 36"H. (Ships KD; FOB Wisconsin, 54174)

ITEM #	SINK DIMENSIONS	TOP
1216K-L*	14"W x 10"D x 6"H	Solid Epoxy

8' INSTRUCTOR'S DESK—Constructed of solid oak, hardwood and oak veneers, with a chemical resistant, earth-friendly UV finish. Included are drawers, and a cupboard for storage of large items. The drawers are expertly crafted from the finest materials with dovetail joinery throughout. Drawers run on epoxy-coated steel glides with nylon bearings for smooth operation. Equipment includes: a hot and cold mixing faucet, double gas cock, GFI protected AC duplex receptacle, polypropylene trap, set of support rods and crossbar, rubber base molding and locks for cupboard and all five drawers (two small, three large). (Ships KD; FOB Wisconsin, 54174)




Small Drawers: 14-3/4"W x 19"D x 5-1/2"H
 Large Drawers: 33-3/4"W x 19"D x 5-1/2"H
 Cupboard: 34-1/2"W x 28"D x 29-1/4"H
 Sink: 18"W x 15"D x 11"H
 Knee Space: 22"W x 32-1/2"H
 Overall: 96"W x 30"D x 36"H

ITEM #	SINK DIMENSIONS	TOP
1116K*	18"W x 15"D x 11"H	Solid Epoxy


Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



ACCESSORIES

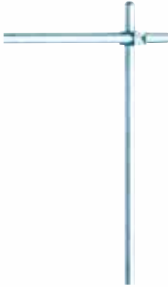

Greenlaw Arm

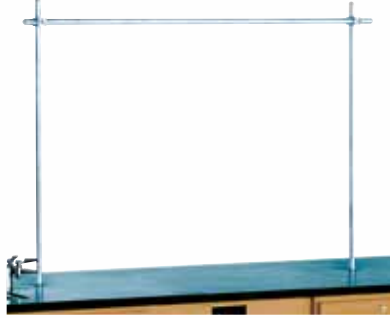

Upright Set

GREENLAW ARM—Includes 3/4" x 3/4" clamp, a 3/4" x 36" upright and a 3/4" x 18" crossbar. Fits in a rod socket with a 3/4" receiver.

UPRIGHT SET—Includes two each 3/4" x 1/2" clamp, two 3/4" x 36" upright and a 1/2" x 40" crossbar. Fits in a rod socket with a 3/4" receiver.

ITEM #	PRODUCT
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar

ITEM #	SIZE	PRODUCT
100002*	3/4" x 3/4"	Clamp
100003*	3/4" x 1/2"	Clamp
100004*	3/4" x 36"	Upright
100005*	3/4" x 18"	Crossbar
100006*	3/4" x 48"	Crossbar
100007*	1/2" x 40"	Crossbar

UPRIGHTS, CROSSBARS, CLAMPS—Solid machined aluminum, these units are nearly impossible to break.

THE LABHAND—Designed to provide a safe, convenient surface for laptop computers in a laboratory setting. The LabHand is mounted on a solid plastic rod which fits securely into standard rod sockets. The platform is made of silver vein colored steel to create a surface that will withstand the typical chemicals found in today's laboratories. It holds laptops and books off the work surface protecting them from spills and increasing the amount of workspace available. Dimensions: 12" square and 5-1/2" rod.



ITEM #	PRODUCT
100994*	LabHand
100994X6*	LabHand (Set of 6)

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

ACCESSORIES



Rod Sockets



Gas Cocks



Gooseneck Faucets



Hand Pump



2-1/2 Gallon Jug



Rubber Leg Boots

ROD SOCKETS—Receives the tapered end of 3/4" uprights.

GAS COCKS—These fittings have a ground key valve for long life. Units are complete with nipple, lock nut and washer for installation and are deck mounted.

GOOSENECK FAUCETS—Chrome plated brass, deck mounted faucets are 9-1/2" H. Multiservice fixtures includes water faucet and two gas cocks. Mounting hardware provided for all units.

HAND PUMP—It has a piston action mechanism for high water flow and long term prime lock. Chrome base and spout.

2-1/2 GALLON JUG—Plastic jug found in our mobile sink units.

RUBBER LEG BOOTS—Fit snugly to 2-1/4" table legs. Sold in set of 4.

FLOOR MOUNTING BRACKETS—Set of four L-shaped mounting brackets are used to permanently attach tables to the floor.

ITEM #	PRODUCT
206505 *	Rod Socket
100047*	Single Gas Cock
100048*	90° Gas Cock
100049*	180° Gas Cock
100466*	360° Gas Cock
100074*	Multiservice Fixture
100051*	Cold Water Faucet
100061*	Hot Water Faucet
100467*	Mixing Faucet
100013*	Chrome Hand Pump
100118*	2-1/2 Gallon Jug
100076X4*	Leg Boots (4)
100728X4*	Floor Mounting Brackets (4)

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



MOBILE DEMO UNITS

EXTRA LARGE MOBILE DEMO CENTER—

Constructed of quality solid hardwood and oak veneers with a chemical resistant, earth-friendly UV finish. It includes one file drawer with full extension glides, 3 smaller drawers (all dovetailed) and two cupboard spaces all with locks and wire pulls. The drawers run smoothly on epoxy-coated steel glides with nylon rollers. It rolls easily on 4" locking swivel casters with a mobile capacity of 500 lbs. Includes: stainless steel sink, hand water pump, two 2-1/2 gallon water bottles and connecting hose, two rod sockets, two each 3/4" diameter x 36"L upright anodized aluminum rods, and 1 each 1/2" diameter anodized aluminum crossbar with connecting clamp, removable storage tote tray and a pegboard with hooks. Also included is the oak framed combination mirror/marker board which has a shatter-proof plastic mirror and can be reversed to serve as a handy marker board. The combination mirror/marker board can be conveniently stored in one of the cupboard spaces. Electrical service is supplied by means of a 25' extension cord which can be connected to a GFI protected AC duplex receptacle installed in the side of the cabinet. Overall dimensions: 30"D x 54"W x 36"H. Shipping weight 290 lbs. (FOB Wisconsin, 54174)



ITEM #	SINK DIMENSIONS	TOP
4342K*	15"W x 13"D x 5"H	1-1/4" ChemGuard



ITEM #	TOP
4501K*	P. Lam
4502K*	ChemGuard

MOBILE DEMO CART—Designed for the transportation of equipment and supplies. Constructed of solid oak with a chemical resistant, earth-friendly UV finish. The legs are 2-1/4" square and are securely fastened to the 4-1/2" high apron with factory installed 5/16" hanger bolts through the heavy-duty 13-gauge steel corner braces. Includes leg stretchers and auxiliary oak plywood shelf. The unit rolls easily on 4 factory installed 3" ball-bearing swivel casters with brakes. The unit is supplied with either a factory installed black 1-1/4" Plastic Laminate or 1-1/4" chemical resistant ChemGuard top. Shipped partially unassembled. Dimensions: 36"W x 24"D x 30"H. Wt. 95 lbs. (FOB Wisconsin, 54174)

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

MOBILE DEMO UNITS



MOBILE LABORATORY UNIT—Features a 1-1/4" thick high-pressure laminated surface with durable drop leaf extensions. Constructed from oak and oak veneers with a chemical resistant, earth-friendly UV finish. It has unbreakable aluminum sliding door glides and locks. The unit comes standard with a stainless steel sink, and hand pump. A GFI protected AC duplex electrical outlet is built into the lab. The cabinet rides on 4" heavy-duty ball-bearing swivel casters with brakes, with a capacity of 500 lbs. Includes a 25-foot extension cord, two tote trays, two 2-1/2 gallon water bottles, connecting hoses, two 3/4" diameter upright rods, one 1/2" diameter crossbar and clamps, a propane burner with stand and 2 adjustable shelves. Weight: 230 lbs. Overall dimensions: 48"W x 24"D x 36"H. Extended Top: 72"W x 24"D. Sink: 15"W x 13"L x 5"D. (FOB Wisconsin, 54174)

ITEM #	PRODUCT
4121K*	Mobile Laboratory Unit
4000K*	Mirror
4001K*	Mirror/Marker Board

MOBILE DEMO TABLE WITH DRAWERS—

A self-contained unit with its own water supply, waste container, stainless steel sink and hand pump, two 3/4" diameter upright rods, one 1/2" diameter crossbar and clamps. The four drawers are mounted on epoxy-coated nylon bearing glides. Sliding, locking doors conceal the drawers and storage compartment for large equipment and the supplied two 2-1/2 gallon water bottles. It is finished with a chemical resistant, earth-friendly UV finish. Heavy-duty 4" ball-bearing swivel casters with brakes, support 500 lbs. A GFI protected AC duplex electrical outlet is built in and can be powered with the supplied 25' extension cord. The top is 1-1/4" thick ChemGuard high pressure laminate. Sink measures 15"W x 13"L x 5"D. Overall dimensions are 48"W x 28"D x 36"H. Weight: 215 lbs. (FOB Wisconsin, 54174)



ITEM #	PRODUCT
4222K*	Mobile Demo Unit
4000K*	Mirror
4001K*	Mirror/Marker Board

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



MOBILE DEMO UNITS

ECONOMY MOBILE LAB TABLE—

Constructed from select oak and oak veneers with a chemical resistant, earth-friendly UV finish. It features unbreakable aluminum sliding door glides and lock. Edges of the 1-1/4" thick high-pressure laminated top are bound with a Plastic-T Molding. With a mobile capacity of 500 lbs., it is equipped with a stainless steel sink, hand pump, heavy-duty 4" locking swivel casters, 25' grounded extension cord, two 2-1/2 gallon water bottles, connecting hoses, two 3/4" diameter upright rods, one 1/2" diameter crossbar with clamps, two fully adjustable shelves, and a GFI protected AC duplex electrical outlet. Sink measures 15"W x 13"L x 5"D. Dimensions: 48"W x 24"D x 36"H. Weight: 185 lbs. (FOB Wisconsin, 54174)



ITEM #	PRODUCT
4111K*	Economy Mobile Lab
4000K*	Mirror
4001K*	Mirror/Marker Board



ITEM #	PRODUCT
4352K*	ADA Mobile Lab Unit
4000K*	Mirror
4001K*	Mirror/Marker Board

ADA MOBILE LAB UNIT—Constructed of solid oak, hardwood and oak veneers with a chemical resistant, earth-friendly UV finish. It is designed as a work station for the physically impaired or for classroom demos with accessibility from both sides. The 1-1/4" ChemGuard top is resistant to most lab chemicals. It includes: stainless steel sink, toggle switch water pump, two 2-1/2 gallon supply and waste containers, GFI AC duplex receptacle with inlet socket, 25' ext. cord, two rod sockets, two each 3/4" diameter by 36" long upright rods, 1/2"x 40" crossbar and two rod clamps. It has a capacity of 500 lbs. and rolls on 4" ball-bearing swivel caster. Dovetailed drawers are constructed of plywood and ride on epoxy coated glides. Cupboard door and drawers have locks. Sink measures: 15"L x 13"W x 5"D. Inside drawers measure: 14.5"W x 19"D x 6"H. Dimensions: 72"L x 24"D x 33"H. Wheelchair access: 32"W x 29"H. Weight: 240 lbs. (FOB Wisconsin, 54174)

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

MOBILE STORAGE



MOBILE DESK WITH STORAGE—The three drawers (one file, two box) have dovetail joints. Features: pull-out writing board, all drawers and doors have locks, 21" x 28" pegboard and hooks for storage, stainless steel sink (15" x 13" x 5"), hand water pump, two 2-1/2 gallon water bottles with all connecting hoses, two 3/4" diameter upright rods, one 1/2" diameter crossbar, two connecting clamps and two rod sockets installed in the table top. A GFI protected AC duplex electrical outlet is built into the lab. Unit has four heavy-duty 4" locking swivel casters. There is a storage area for the optional mirror/marker board. It is finished with a chemical resistant, earth-friendly UV finish. The top is a 1-1/4" chemical resistant ChemGuard top. Weight: 250 lbs. (FOB Wisconsin, 54174)

ITEM #	PRODUCT	SIZE
4332K*	Mobile Desk w/Storage	28"D x 48"W x 36"H
4000K*	Mirror	34-1/2"W X 22-1/2"H
4001K*	Mirror/Marker Board	27-3/4"W x 20-3/4"H

MOBILE STORAGE CABINET—A compact, rolling storage cabinet that is constructed of select oak and oak veneers with a chemical resistant, earth-friendly UV finish. The adjustable shelf allows tall items to be placed in the unit as well as allowing use of the shelf (measures: 32"W x 18"D). With a mobile capacity of 500 lbs., it comes equipped with four heavy-duty swivel casters with brakes. The doors lock for security. The unit is available with black 1-1/4" plastic laminate or 1-1/4" ChemGuard top. Shipping weight: 140 lbs. (FOB Wisconsin, 54174)



ITEM #	TOP	SIZE
4401K*	P. Lam	36"W x 24"D x 36"H
4402K*	ChemGuard	36"W x 24"D x 36"H

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



MOBILE STORAGE

BARE BONES MOBILE STORAGE CABINET—

Constructed of select oak and oak veneers with a chemical resistant, earth-friendly UV finish. Includes unbreakable aluminum sliding door glides, lock and one fixed shelf with a clearance of 13-10/16. Edges of the 1-1/4" thick high pressure laminated top are bound with a Plastic-T Molding. With a weight capacity of 500 lbs., it rolls on 4" heavy-duty ball-bearing swivel casters with brakes. Weight: 180 lbs. (FOB Wisconsin, 54174)



ITEM #	PRODUCT	SIZE
4101K*	Bare Bones Mobile Storage Cabinet	48"W x 24"D x 36"H



MOBILE TOTE TRAY STORAGE—

Constructed of oak veneers and solid oak with a chemical resistant, earth-friendly UV finish. With a mobile capacity of 500 lbs., this unit rolls easily on 4" locking swivel casters. The two lockable doors swing 270° and conceal 21 tote trays (19"L x 14-1/4"W x 3-1/2"H). Plastic laminate top can double as a work surface. Overall unit measures: 48"W x 24"D x 41-1/2"H. (FOB Wisconsin, 54174)

ITEM #	PRODUCT
4751K*	Mobile Tote Tray Storage
100135*	Replacement Tote Tray

MOBILE BALANCE STORAGE CABINET—

Constructed of solid oak and oak veneered plywood with chemical resistant, earth-friendly UV finish. Featured are locking doors, both front and back with unbreakable aluminum glides. Supplied with four adjustable shelves (weight capacity of 40 lbs per sq. ft.). With a mobile capacity of 500 lbs., it rolls on heavy-duty swivel casters with brakes. The top is black 1-1/4" plastic laminate and furnishes an auxiliary work surface. Dimensions; Shelves are 3/4"H x 8"D x 42"L, Overall is 48"W x 24"D x 40"H. Weight: 230 lbs. (FOB Wisconsin, 54174)



ITEM #	PRODUCT
5201K*	Mobile Balance Cabinet Storage

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

MOBILE STORAGE



MOBILE MICRO-CHARGE STATION—

A mobile oak cabinet that is equipped to recharge any **Ken-A-Vision** cordless microscope within eight hours while it is safely stored behind locked doors. The peninsula style cabinet interior (7"W x 9"D x 15"H) holds 16 microscopes and has storage space for other essentials such as slides. The LED meter lets you know when the microscopes are charging and when they are fully charged. All this with only one electrical connection from the cabinet's side. Is is finished with a chemical resistant, earth-friendly UV finish. The 1-1/4" plastic laminate top doubles as a work surface. Mobile capacity: 500 lbs. Dimensions: 48"W x 30"D x 41-1/2"H. Weight: 311 lbs. (FOB Wisconsin, 54174)



ITEM #	PRODUCT
4741K*	Mobile Micro-Charge Station



MOBILE MICROSCOPE STORAGE CABINET—

Provides a safe, convenient way to store your microscopes. It is constructed of solid hardwood, solid oak and oak veneers with a chemical resistant, earth-friendly UV finish. The front and back doors have locks. The unit contains removable wire dividers to provide an individual space for each

individual microscope (6" x 11" x 16"). With a mobile capacity of 500 lbs., this cabinet rolls on heavy-duty swivel casters with brakes. The top is 1-1/4" black plastic laminate and furnishes an auxiliary work surface. The cabinet has storage for 24 standard classroom microscopes. Weight: 220 lbs. (FOB Wisconsin, 54174)

ITEM #	PRODUCT	SIZE
4701K*	Mobile Microscope Storage Cabinet	48"W x 24"D x 40"H


Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



TALL STORAGE

TALL STORAGE W/GLASS DOORS—

Showcase your specimens with this beautiful clear glass door cabinet. Solid hardwood, solid oak and oak veneers, protected by a chemical resistant, earth-friendly UV finish. The oak framed 3/16" tempered glass doors are designed to crumble when broken preventing large dangerous shards. Each cabinet comes with a three-point locking system. Glass door units have two fixed and four adjustable 1" shelves (weight capacity: 40 lbs. per sq. ft). (FOB Wisconsin, 54174)

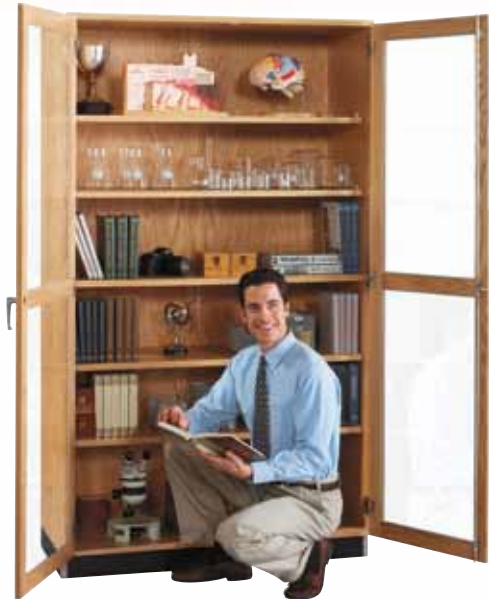

ITEM #	PRODUCT	SIZE
358-4822*	48"W Glass Door Tall Storage	22"D x 84"H

Cabinet should be secured to the wall for safety. Mounting hardware is not included.



TALL STORAGE W/OAK VENEER DOORS—

Solid hardwood, solid oak and oak veneers with a chemical resistant, earth-friendly UV finish. Available with hinged oak doors. Each cabinet comes with a three-point locking system. All units include two fixed and four adjustable shelves. Weight capacity: 40 lbs per sq. ft. (FOB Wisconsin, 54174)

ITEM #	PRODUCT	SIZE
353-4822*	48"W Oak Door Tall Storage	22"D x 84"H

Cabinet should be secured to the wall for safety. Mounting hardware is not included.

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

TALL STORAGE



TOTE TRAY STORAGE CABINET—

This organized storage case holds 48 tote trays (tray dimensions:) which are included. A great place for students to store their art or science projects. Tote tray storage case has hinged oak doors with three-point locking handles. It is finished with a chemical resistant, earth-friendly UV finish. Weight: 450 lbs. (FOB Wisconsin, 54174)

ITEM #	PRODUCT	SIZE
350-4822*	Tote Tray Storage Cabinet	48"W x 22"D x 84"H
100135*	Replacement Tote Tray	19"W x 14-1/2"D x 3-1/2"H

Cabinet should be secured to the wall for safety. Mounting hardware is not included.

MICROSCOPE STORAGE CABINET—Constructed of solid oak, hardwood and oak veneers, using dowel construction, with a chemical resistant, earth-friendly UV finish. Solid oak framed glazed doors are hinged with three-point locking handles. Tempered glass is 3/16" thick. Microscopes are stored safely and conveniently on two fixed and three adjustable 1" shelves (weight capacity of 40 lbs per sq. ft.). Adjustable wire dividers allow for storage of microscopes of various sizes. Dividers can be removed to use shelves for general storage. Holds 30 microscopes, individual compartment measures: 7-1/4"W x 14"D x 14-1/2"H. (FOB Wisconsin, 54174)



ITEM #	SIZE
372-4816*	48" W x 16"D x 84"H

Cabinet should be secured to the wall for safety. Mounting hardware is not included.

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



PLANNING TABLES



Maple



Almond P. Lam



Natural Birch P. Lam

PLANNING TABLES—Our multi-use table with a base constructed of two tapered solid maple pedestal style panel legs and two solid maple truss bars. The table is finished with an earth-friendly UV finish. Choose from 1-3/4" solid maple top or 1-1/4" plastic laminate with 3 mm PVC edging in either birch or almond top colors. All units are 60"W x 42"D.

ITEM #	SIZE	TOP
PT-62M*	60"W x 42"D x 30"H	Maple
PT-62PNB*	60"W x 42"D x 30"H	Natural Birch P. Lam
PT-62P*	60"W x 42"D x 30"H	Almond P. Lam



THE WORKTOP CLASSIC—

Designed to provide a structurally sound and completely functional work center. This table is constructed of solid maple legs and apron with 13-gauge steel corner braces. The unit has solid 2-1/4" maple legs which are outfitted with adjustable floor glides. The table is finished with an earth-friendly UV finish. The table comes standard (30"H) height and with an almond colored 1-1/4" plastic laminate or a 1-3/4" maple top.

ITEM #	SIZE	TOP
WX4-M*	48"W x 48"D x 30"H	Maple
WX4-P*	48"W x 48"D x 30"H	Almond P. Lam

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

PLANNING TABLES



SHAIN DRAWING TABLE—Your choice of 12 different tables includes one- and two-piece tops, plain apron, two sizes of drawers and two heights. All units have solid 2-1/4" maple legs and aprons with fully adjustable 3/4" almond colored plastic laminate tops with the soft close feature. The one-piece tops measure 36"W x 24"D. The two-piece top has a 24"W x 24"D adjustable side and a 12"W x 24"D fixed, flat side. The DT-1 models have a large center drawer (20-1/2"W x 18"D x 3"H). The DT-2 and DT-3 models have a small drawer (9"W x 18"D x 3"H) on the right side of the apron. The DT-9 model features a plain apron. The units are finished with an earth-friendly UV finish. All models are available in either 30"H or 36"H.

APRON MODEL	1-PIECE TOP (30")	2-PIECE TOP (30")	1-PIECE TOP (36")	2-PIECE TOP (36")
Plain	DT-9A30*	DT-9SA30*	DT-9A37*	DT-9SA37*
Small Drawer	DT-2A30*	DT-2SA30*	DT-3A37*	DT-3SA37*
Large Drawer	DT-1A30*	DT-1SA30*	DT-1A37*	DT-1SA37*

ART BENCH—Made of 3/4" solid maple. A built-in shelf provides storage. Finished with an earth-friendly UV finish.



ITEM #	PRODUCT	SIZE
AB*	17" Bench	17"W x 12"D x 17"H



ADAPTABLE TABLE (ADA)—The height can be adjusted easily. The unique design of the leg allows you to adjust the height of the table in less than 30 seconds without the use of tools. Use the table at sitting, standing, or display height. Adjust only one side of the table to create an easel. Ideal for use by a person in a wheelchair. Make it match the height they need. Table adjusts in 2 cm increments. The powder-coat legs include adjustable glides. Choose from two sizes. Height adjusts from 26-39"H.

ITEM #	TOP SIZE	TOP
ALT-6030BL*	60" x 30"	Black P. Lam
ALT-6042BL*	60" x 42"	Black P. Lam

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet



DRAFTING TABLES

DRAFTING/DRAWING TABLE—Our top selling table for drafting and computer generated art use. Designed for rooms with limited space options, this table is perfect for drafting, graphic art and CAD. The solid maple pedestal legs support this sturdy table. A rear shelf (9-7/8"D x 42"W) can hold a computer monitor. The 42"W x 30"D adjustable top has a soft close feature and includes a pencil stop. Options include: Flat panel monitor arm, stationary monitor platform, and factory installed under-desk keyboard tray. All units are 30"H.



CDTC-1VT M UDK

TOWER MODEL	OPTIONS
CDTC-1VT*	Table Only
CDTC-1VT UDK*	w/Keyboard Tray
CDTC-1VT SP*	w/Stationary Platform
CDTC-1VT M*	w/Flat Panel Monitor Arm
CDTC-1VT SP UDK*	w/Keyboard Tray & Stationary Platform
CDTC-1VT M UDK*	w/Keyboard Tray & Flat Panel Monitor Arm



DT-60SA



CDTC-70

DRAWING TABLE—Constructed of solid maple framing and reinforced with 13-gauge steel corner braces. Your choice of a one- or two-piece top made of 3/4" almond colored plastic laminate and equipped with a soft-close feature and a pencil stop. On the one-piece top models, the fully adjustable top measures 42"W x 30." On the two-piece top, the 30" x 30" left side of the top can be fully adjusted and is equipped with a soft-close feature and a pencil stop and the 12"W x 30"D flat top on the right side. The deluxe model includes a CPU holder, a monitor arm, keyboard & mouse tray, and cable manager. Overall size: 42"W x 30"D x 39-3/4"H.

ITEM #	TOP STYLE	TOP
DT-30A*	Basic, One-Piece	Plastic Laminate
DT-60SA*	Basic, Two Piece	Plastic Laminate
CDTC-70*	Deluxe, Two-Piece	Plastic Laminate

Quantity Restrictions Apply: 30 Per Table, 6 Per Cabinet

WORKBENCHES



OPEN-STYLE WOOD WORKBENCH—

It has a large work surface and open storage space on the maple plywood lower shelf. Available with or without vises. The base is constructed of solid 2-1/4" maple legs and frame. The top is 2-1/4" thick maple; the top size is 54" x 64".

ITEM #	# VISES
WW4-0V*	0
WW4-4V*	4

STEEL WORKBENCH—Base consists of sturdy heavy gauge steel legs with stringers. Choose between two-station and four-station units. Top is constructed of 2-1/4" maple. Two-station measures: 64"W x 28"D x 32-1/4"H. Four-station measures: 64"W x 54"D x 32-1/4"H.



ITEM #	STYLE
WBML2-0V*	Two-Station
WBML2-1V*	Two-Station w/1 Vise
WBML2-2V*	Two-Station w/2 Vises
WBML4-0V*	Four-Station
WBML4-4V*	Four-Station w/4 Vises

FOUR-STATION STEEL WORKBENCH—

Base consists of a double-faced heavy gauge steel unit, welded and riveted throughout. The unit provides 12 locker openings (12"W x 21"D x 15"H). Vented doors have spring loaded hinges, padlock hasp, and knockout plug for cylinder locks. Top is constructed of 2-1/4" maple. Overall size: 64"W x 54"D x 33-1/4"H.



ITEM #	STYLE
WB12-0V*	Gray
WB12-4V*	Gray w/4 Vises